


it's time to change your point of view

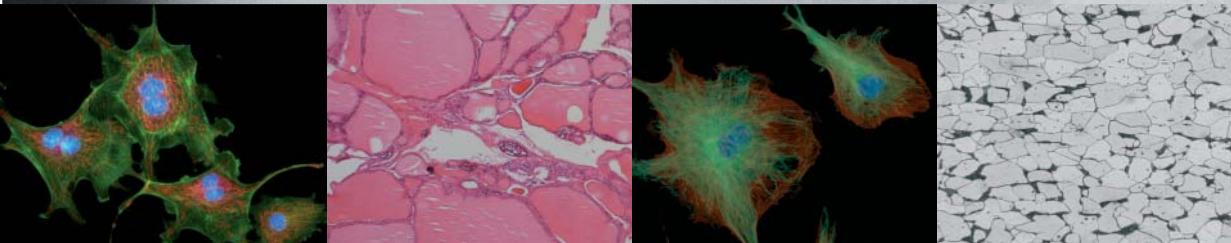

ProgRes® 14bit cameras Monochrome and Color

With the new monochrome camera series **ProgRes® MF** JENOPTIK offers a professional solution for the special tasks in fluorescence microscopy.

The 14bit cameras with a resolution of 1.4 million pixel deliver detailed images and assure a sparing treatment of your fluorescence samples due to high image refresh rates of up to 51 frames per second offered in several binning-modes. The sensitivity and the quality of the image is additionally improved by an analogue gain and a switchable thermo-electrical cooling, optionally available.

A sophisticated fluorescence mode, additionally integrated

in the **ProgRes® Capture Software**, offers multi-color fluorescence image capturing and enables the user to minimize auto-fluorescence already during image capture.

Additionally available as color camera series **ProgRes® CF**, the systems are also suitable for applications in quality assurance, as the 14bit digitization delivers superb contrasts and high information depth required for material analysis.

Of course, the cameras can be externally triggered and are easily adoptable to any microscope or PC in combination with C-Mount and FireWire® interface.

Camera Models:

ProgRes® MF

- *Monochrome CCD camera, 12bit, non-cooled*

ProgRes® MF^{cool}

- *Monochrome CCD camera, 14bit, cooled*

ProgRes® CF


- *Color CCD camera, 12bit, non-cooled*

ProgRes® CF^{cool}

- *Color CCD camera, 14bit, cooled*




ProgRes®

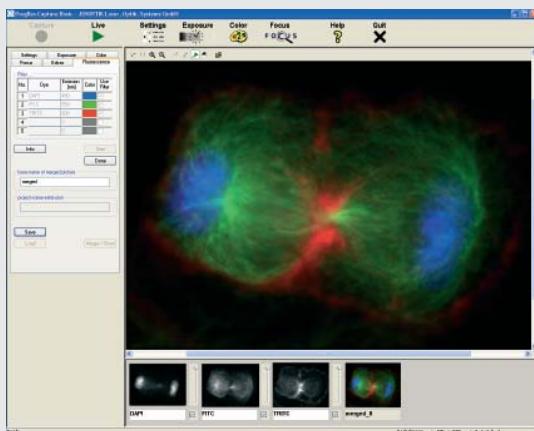

Fluorescence mode for capture of multi-color florescences – user-guided capture procedure and combination of several filter images.

Application Areas:

- Phase contrast
- Fluorescence
- Transmitted & Reflected Light
- Brightfield & Darkfield

In the fields of:

Live-Science

- Genetics
- Microbiology
- Pathology
- Histology
- Hematology
- Cytology
- Zoology

Material sciences

- Quality control
- Geology
- Mineralogy
- Chemistry
- Semiconductor industry
- Industrial inspection

Forensics

Technical Data

Camera Model	ProgRes® MF	ProgRes® MF ^{cool}	ProgRes® CF	ProgRes® CF ^{cool}
CCD Sensor	2/3" 1.4 Megapixel CCD Sensor (Sony ICX285AL/AQ) Active area: 8.8 x 6.6 mm ²			
Cooling	no	Peltier & Fan	no	Peltier & Fan
Digital output	12bit	14bit	12bit	14bit
Sensor resolution	1360 x 1024			
Pixel size	6.45 µm x 6.45 µm			
Read-out frequency	Switchable: 12 MHz and 24.5 MHz			
IR cut-off filter	Hoya C-500S (optionally available without)			
Dynamic range	ca. 66 db/typical: >2000:1			
Exposure time	0 to 600 s			
Gain	Analogue 1x to 8x			
Image refresh rate	Up to 51 fps			
Image Resolution	Programmable Resolution 1360 x 1024 680 x 512 High Frame Read-out Mode 340 x 256 High Frame Read-out Mode 680 x 512 2x2 binning (Fluorescence) 452 x 340 3x3 binning (Fluorescence) 340 x 256 4x4 binning (Fluorescence) 272 x 204 5x5 binning (Fluorescence) FireWire® IEEE1394a C-Mount (0.63x TV Adaptor for microscope usage) External trigger in and out Dual thread 3/8" and 1/4"			
Digital interface				
Optical interface				
Trigger				
Tripod thread				
Software	ProgRes® Capture Software for MS Windows® 2000/XP (TWAIN and Stand-Alone) ProgRes® Capture Software for Apple Macintosh OS X (in preparation) SDK for MS Windows® and Apple Macintosh OS X PC: Pentium IV 1.4 GHz or better; ≥512 MB RAM; FireWire® (OHCI Standard) MAC: G4 or better; 512 MB RAM			
Hardware requirements				
Power consumption	6 W	8 W	6 W	8 W
Weight	800 g / 1.76 lbs			
Dimension	145 x 93 x 123 mm / 5.7 x 3.6 x 4.85 " (L x W x H)			
Operating conditions	Temperature: +5°C to +35°C (41°F to 95°F) Humidity: 5% - 80%, not condensing			

This design and related specifications are subject to continuously ongoing development. We reserve the right to make changes in the interest of technical progress.

JENOPTIK
Laser, Optik, Systeme GmbH
Business Unit Sensor Systems
07745 Jena, Germany
Phone +49 3641 65 39 63
progres@jenoptik.com
www.progres-camera.com

Your direct sales agent for high-grade microscope cameras: